


HTM02-01 Chapter 3

Mounting height range for AVSU's

1000 - 1800mm (1000mm is the optimum height: when multiple-circuits are installed, it may be necessary to reduce this to 800mm to avoid too many columns.)

Mounting order

- O2
- N2O
- O2/N2O
- MA4
- SA7
- VAC
- AGSS

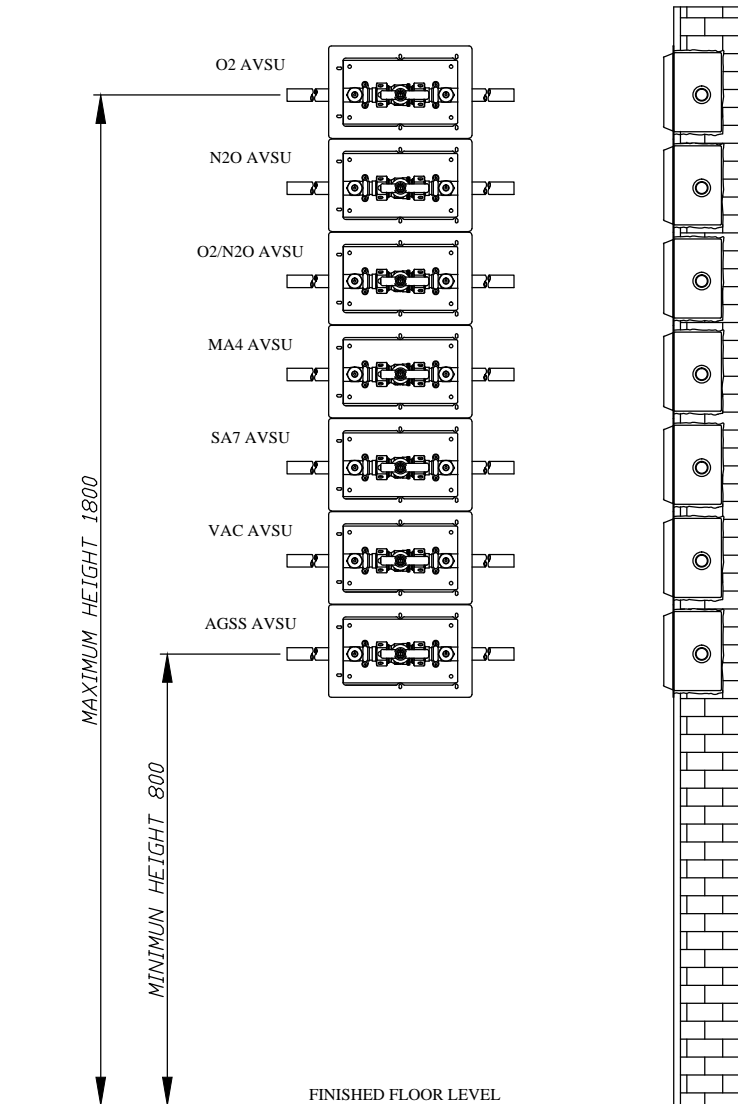

If In Doubt Ask							
Do Not Scale	Drawn by K.J.P.	Checked by	Approved by – date	Filename	Date 21/01/2014	Scale NTS	All Dims In mm Unless Otherwise Stated
<p><b>PRECISION UK Ltd</b>                  CPX Technology Building, Pepper Road, Hazel Grove, Sockport SK7 5BW                  Tel: 00 44 (0) 161 487 2822 Fax: 00 44 (0) 161 487 2816                  Email: info@precisionuk.co.uk www.precisionuk.co.uk</p>				AVSU INSTALLATION DETAIL			
				AVSU-INST-02	Revision 0	Sheet 1/1	